


SendToPACS

CONVERT NON-DICOM FILES TO DICOM AND SEND TO PACS



VERY EASY TO USE

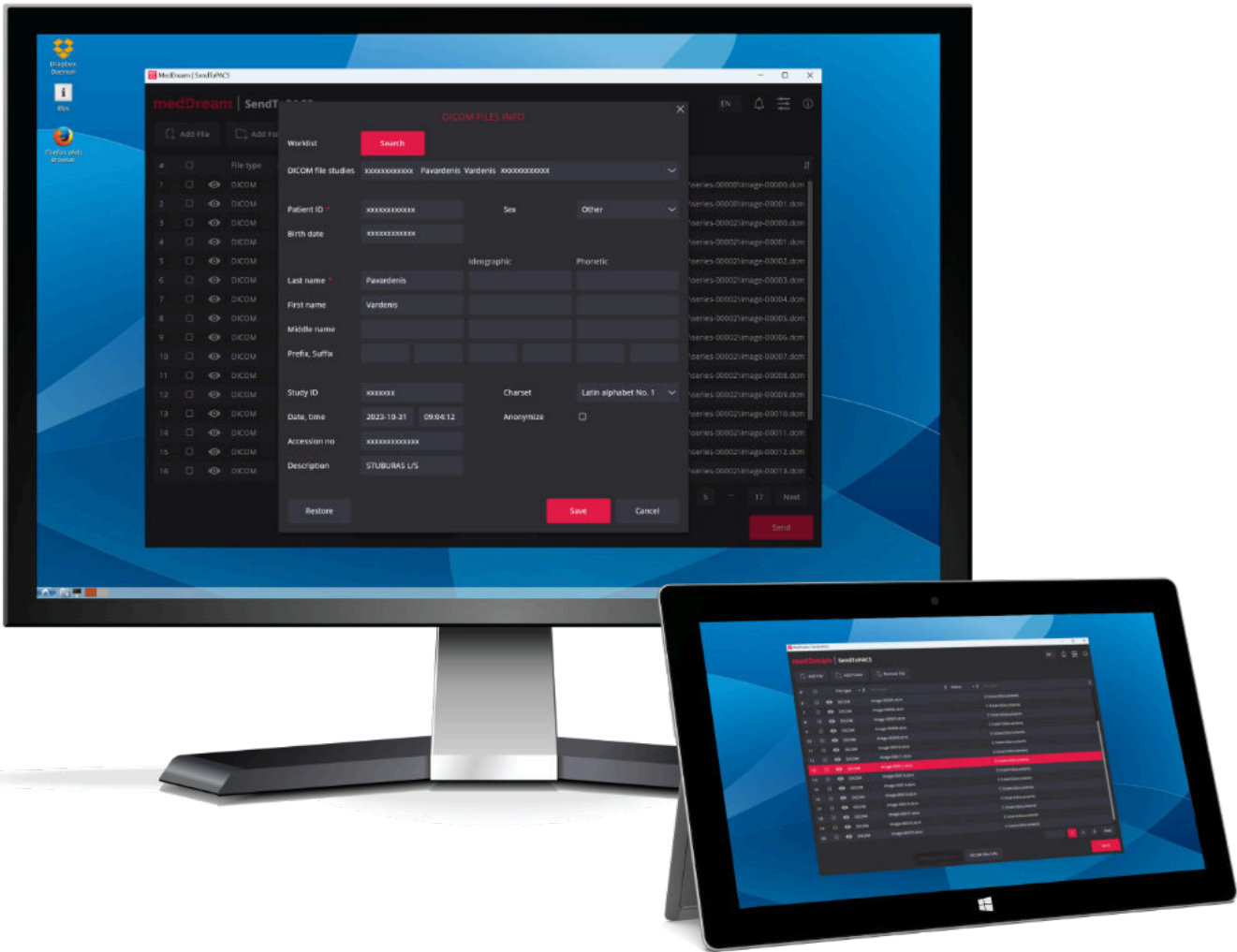

IMAGES, VIDEOS, SIGNALS AND DOCUMENTS SUPPORT



ENHANCED CHARACTER SET SUPPORT




DICOM MODALITY WORKLIST


medDream

Non-DICOM files to PACS

SendToPACS software converts non-DICOM images and video files into DICOM format files, and sends them to PACS server or DICOM workstation.

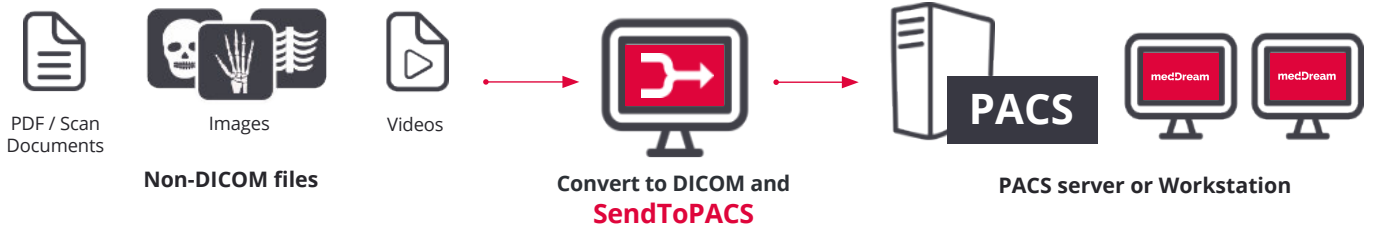

DICOM files to PACS

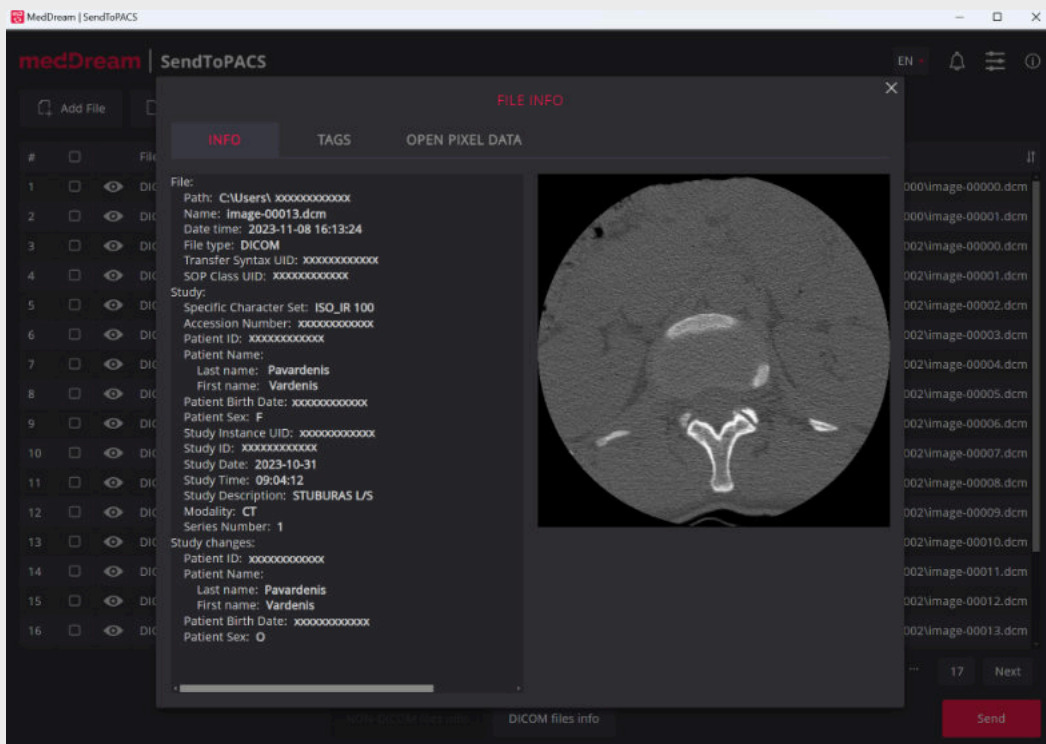

SendToPACS features:

Main features:

- Files conversion to DICOM with the highest possible quality;
- Selectable video conversion to DICOM as MPEG2 or MPEG4 format;
- Sending DICOM files to PACS server or DICOM workstation;
- DICOM modality worklist support, advanced worklist search;
- Connectivity verification (echo);
- Windows context menu support;
- Drag and drop files or folders;
- Multi-language support.

Customization:

- Displays worklist result in UTF-8 or Latin1;
- Possibility to query HTTP/HTTPS worklist;
- Customizable modality representation value;
- DICOM specific character set latin1 (ISO_IR 100) or UTF-8 (ISO_IR 192);
- Customizable SOP Class UID;
- Customizable video conversion process CPU resources (number of threads).

Installation requirements:

- Supports Windows OS.

SendToPACS supports the following files and formats

Images:

jpg, png, tif, bmp, dcm.

Videos:


mov, mpeg, mpg, avi, mp4, mov, wmv.

Documents:



Structured Report and PDF files.

Signals:

ECG and others.

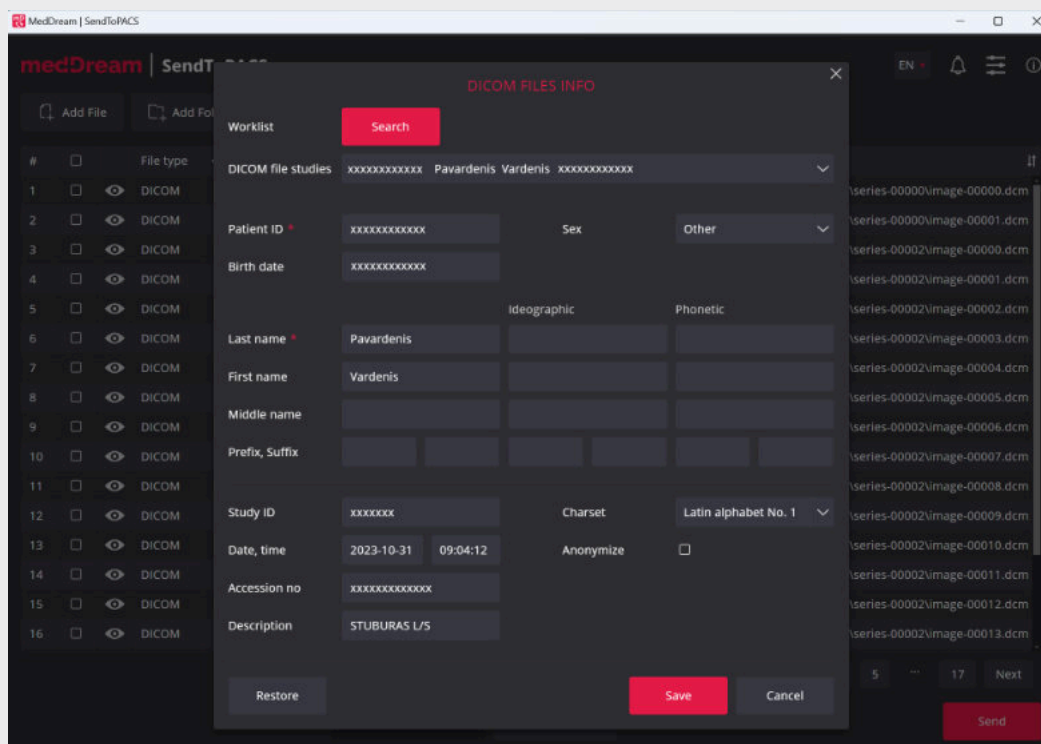

Workflow with SendToPACS

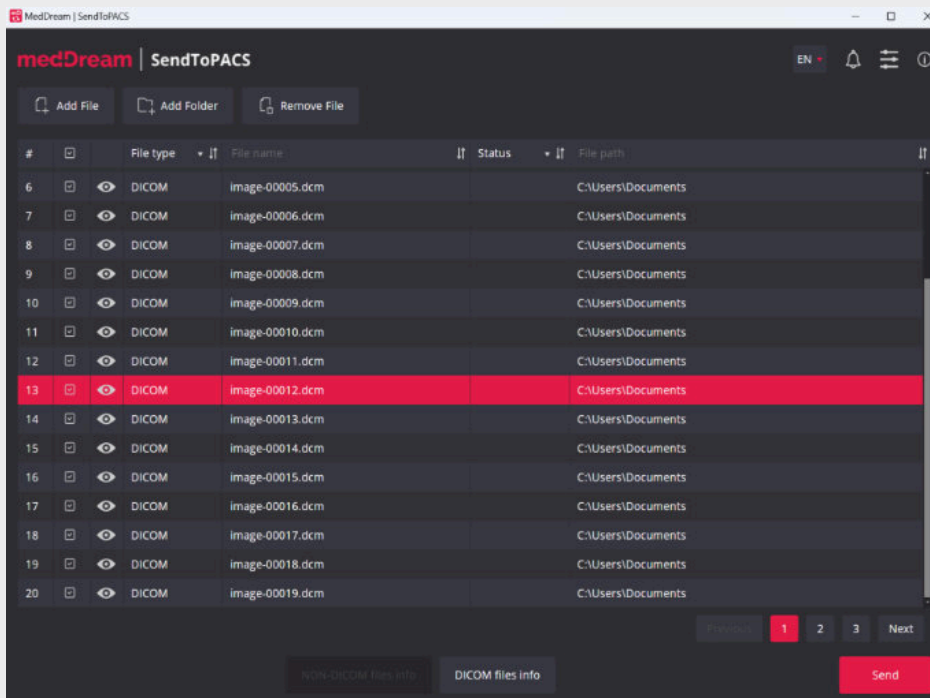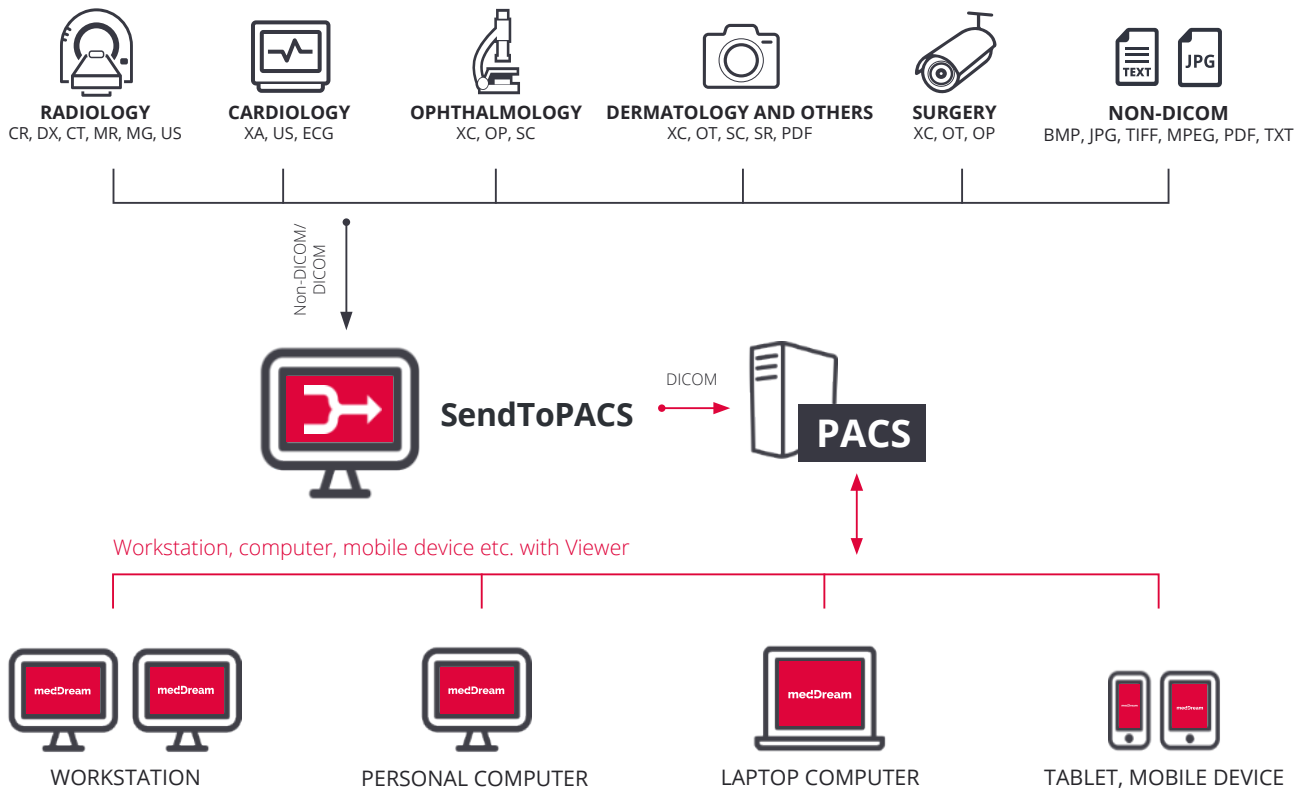

medDream

E-mail: info@meddream.com
 Phone: +370 670 15991
 Web: www.meddream.com

Address: K. Baršausko St. 59,
 LT-51423, Kaunas,
 Lithuania